

# Mini Patient Hand Units

Issue 3 (02/21) - 1



STATIC SYSTEMS  
HEALTHCARE  
COMMUNICATION SOLUTIONS



## Overview

Stylishly contoured to fit comfortably in the hand, the innovative design of Ultima and Codemlon hand units enables patients to quickly identify their needs and make confident choices, whilst the robust construction of the units ensures suitability for use within the healthcare environment.

This product can be used with the following system(s):

- Ultima**
- Codemlon**
- Aspire SmartSync**

## Features

- **Ergonomic design and distinguishing colours** as recommended by the RNIB for ease of patient use
- **Sealed to IP 67** for improved product reliability and to allow disinfecting between patients for improved infection control
- **Tactile buttons** for user reassurance
- **Nurse call button with permanently illuminated back light** for easy locating in the dark
- **Reassurance LED** illuminates to acknowledge operation of the nurse call button
- **LED torch facility** to illuminate the immediate bed area without disturbing adjacent patients
- **Quick-release plug** helps prevent damage to equipment
- **Cable strain relief** to minimise potential cable damage
- **Replaceable cable** extends product life
- **Linen clip** to prevent the unit slipping out of reach
- **High immunity to electro-static discharge** to protect patient, staff and equipment
- **Co-moulding** for improved grip and seal
- **HTM Compliant**

## Part No. Description

- 75801B** Mini Patient Hand Unit - lighting
- 75821B** Mini Patient Hand Unit - lighting & catering
- 75875B** Mini Patient Hand Unit - lighting
- 75876B** Mini Patient Hand Unit

*Note: Hand Units are supplied without a Mounting Bracket (Part No. 37140B).*

### Replacement Parts List:

- 42942 Replacement Cable with linen clip and plug

*Note: These units do not contain any user serviceable parts, however Static Systems offer a service exchange facility and it is possible to replace the cable.*




Static Systems Group Limited  
Heath Mill Road, Wombourne,  
Staffordshire, WV5 8AN, UK  
T: +44 (0) 1902 895 551



A Halma company

## Dimensions

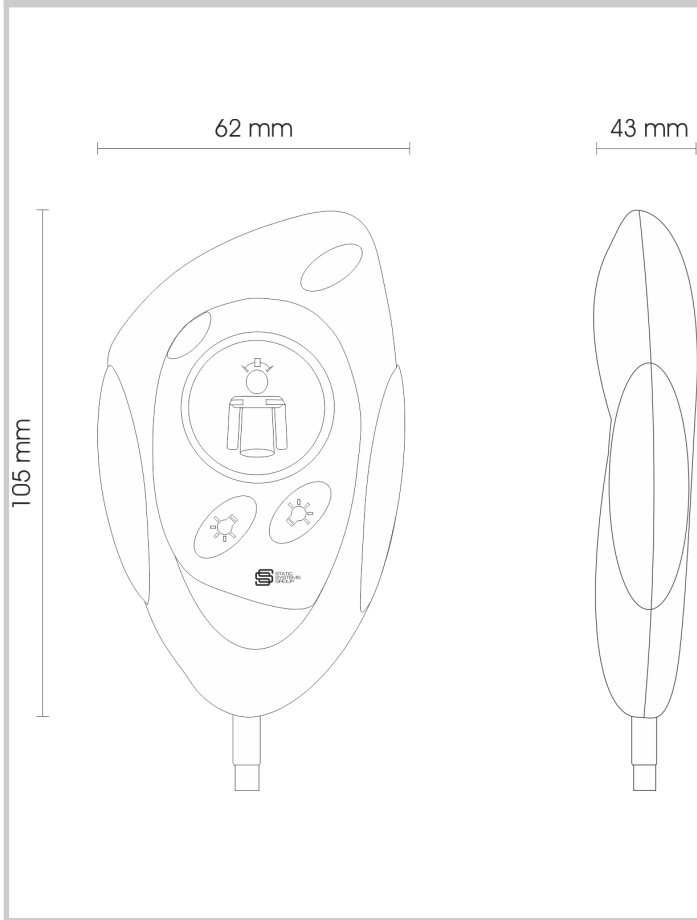

## Specification

Material: Silicone Rubber, ABS, TPE

Weight: 106g

## Features

	75801W	75821W	75875W	75876W
Reassurance LED	✓	✓	✓	✓
Patient Call Push	✓	✓	✓	✓
Catering Call Push		✓		
Lighting Control	✓	✓	single ✓	
LED Torch Facility	✓	✓	✓	✓
Linen Clip	✓	✓	✓	✓
Quick Release Plug	✓	✓	✓	✓
2 Metre Flexible Cable	✓	✓	✓	✓

## Bedlight Control

Where hand units include on/off functionality for the control of bed lighting, the buttons can be programmed to allow variable dimming. Appropriate electronics must also be provided in the lighting unit.

Bed light buttons will only operate when used in conjunction with Static Systems' bedhead services and nurse call panels that comprise a bedlight control section with relays. When hand units are used with a basic nurse call unit such as a reset/jack/emergency panel, the bedlight control function will not be available.